

Technical Data Sheet

Doc No–Ab/TDS/BC-148

Rev No: 05 Date 04/2023-24

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Amber Wide Mouth Bottle Graduated, HDPE

Valid for Catalog No. P11180A, P11181A, P11182A, P11183A, P11184A,P11185A

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

Amber Wide Mouth Bottle , HDPE

- 1.1.1 These bottles come with inner thread for superior sealing and better security.
- 1.1.2 Opaque amber colour reduces UV light transmissions to protect light sensitive samples.
- 1.1.3 Wide mouth of bottles allows for easy filling and emptying.
- 1.1.4 These bottles are used for packing, shipping and storage purposes.
- 1.1.5 Superior chemical resistance.
- 1.1.6 Complete emptying is guaranteed through inclination of shoulders.

1.2 Raw material:

1.2.1 Body: HDPE (High Density poly Ethylene)

1.2.2 Cap: PP (Polypropylene)

1.3 Material:

- 1.3.1 High Density poly Ethylene (HDPE), Hydrophobic, confirm US FDA 21CFR
- 1.3.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for pharmaceutical use.
- 1.3.3 The product has passed the USP testing including class VI tests, Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.
- 1.3.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE ,2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water.
- 1.3.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.4 Quality Control:

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (ml)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11180A	30	21.2	58.9	72
P11181A	60	24.8	81	72
P11182A	125	32	94.5	72
P11183A	250	32.2	126.7	72
P11184A	500	47.5	163.2	48
P11185A	1000	61.5	196.8	24

1.6 GENERAL FEATURES:

CE	Yes
Autoclavability	No
Translucence	Yes
Resistant	To be effected by radiations
Sterilization	Yes
Temperature	From – 40°C to + 80°C
Leak Proof	Yes

1.7 Chemical resistance: See website (www.abdoslifesciences.com)

1.7.1 Manufactured in Aseptic Environment

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.

1.8.3 Poor resistance; exposure to chemical is not allowed or causes immediate permanent damage. Results of testing with reference materials are to be considered as indication: in case of specific use, it is recommended to carry out preliminary testing.

Information relating to the resistances of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products.

1.9 Sterilization:

Sterile catalog number from this lot have been pre-sterilized and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products Requirements for validation and Routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation.

WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure.

Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.

Consult OSHA guidelines prior to handling or lifting filled carboys/bottles/aspirators.

If applicable: Using the spigot as a handle can ruin the threads and causing the spigot not to seal properly.

Cleaning:

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

Cap and Spigot Cleaning: (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA



Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslivesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

UKRP

Obelis UK LTD
Sandford Gate,
East Point Business Park,
Oxford OX4 6LB
United Kingdom
Tel: +44.1491.378012
info@obelis.co.uk

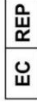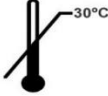

Obelis s.a.
Bd General Wahis 53
1030 Brussels Belgium
Tel: +(32) 2 732-59-54
Fax: +(32) 2 732-60-03
mail@obelis.net



Abdos Labtech Pvt. Ltd.

Plot No. 1 Lakshwari Bhagwanpur Roorkee 247667 (UK) India
Email : labtech@abdosindia.com, www.abdoslifesciences.com
Ph: +(91) 11 41081897



Electronic IFU

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakshwari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903

Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7