

Technical Data Sheet

Doc No – Ab/TDS/PCR-108

Rev No: 07 Date 04/2023-24

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Sealing Mats, Silicon Rubber

Valid for Cat. No. P10721, P10731

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Sealing Mats, Silicon Rubber.

1.1.1 These mats are made of silicone rubber and designed to seal the plates

1.1.2 The mat is made up of 96 plugs & 384 plugs. The mat has to be compressed. A plug that caps the individual well seals each well on the PCR plate.

1.1.3 Sealings mats are autoclavable.

Design:

Sealing Mats designed to seal all PCR plates.

1.3 Performance Specification:

Code	Type	No. per Case
P10721	Silicone Sealing Mats 96 well	50 mats
P10731	Silicone Sealing Mats 384 well	50 mats

1.4 Material:

1.4.1 Silicon Rubber, confirm US FDA 21CFR

1.4.2 Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.3 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE, 2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water

1.4.4 Product is manufactured without the use of Raw Materials of animal origin Therefore BSE/TSE issue does not concern this product to the best of our knowledge.

1.5 Color:

Translucent (Natural)



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email: labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.6 Quality control:

Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.7 Other information

1.7.1 DNase, RNase, DNA, Endotoxin and PCR inhibitor.

1.7.2 Recommended for Single use only

1.7.3 These mats are designed to seal all PCR plates except P10717.

1.8 Basic Features

1.8.1 Temperature range: -80°C to + 120°C

1.8.2 Synthetic elastomer in which silicon atoms are linked together by oxygen atoms.

1.8.3 Very good thermal stability and resistance to UV, Ozone and weather

1.8.4 Autoclavability: Autoclavable.

1.8.5 Manufactured in Aseptic Environment

1.8.6 Chemical resistance: see website (www.abdoslifesciences.com)

*The expiry of Sealing Mats without glue is generally 3 years. Because there is glue on the Mats, and company's storage environment will affect the stickiness of the glue, the default validity period is one year.

Autoclaving Guideline for Plastics

All plastic consumable are Autoclavable at 121°C for 20 minutes, 15 psi/1atm. Do not use dry cycle as this may cause deformation or weakening of plastic. Allow temperature within the autoclave to return to at least 50°C before removing product.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.

WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.




Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakshewari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email: labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

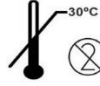

Abdos Labtech Pvt. Ltd.

Plot No. 1 Lakshwari Bhagwanpur Roorkee 247667 (UK) India
Email : labtech@abdosindia.com, www.abdoslifesciences.com
Ph: +(91) 11 41081897




Electronic IFU

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110

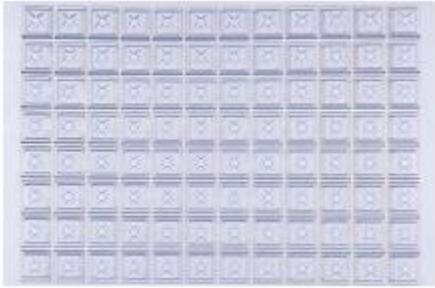

Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakshwari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903

Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7