

TECHNICAL DATA SHEET

Doc No – Ab/TDS/Cryo-27
Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Cryo Vial Internal Threaded Sterile with Silicon Seal, PP

Valid for Cat. No. P60108, P60109, P60110, P60112, P60113

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Disposable Cryo vials, DNase, RNase, Human DNA, ATP & Pyrogen Free, IVD certified

1.2 Design:

1.2.1 Vials: Available in Star Foot, Round Bottom with Internal Thread

1.2.3 Cap: Assembled Cap (PP) along with Silicon Washer

1.3 Performance Specifications:

1.3.1 Vial Volume: 1.0 ml, 1.8 ml and 4.5 ($\pm$ 0.1 ml)

1.3.2 Vial Size:

Code	Capacity (ml)	Outer Diameter (mm)	Height	No. Per Case	Self-Standing	Bottom	Brimful Capacity (Approx ml)
P60108	1.0	12.3	41.2	1000 (25x40)	Yes	Star Foot	1.70
P60109	1.8	12.3	47.5	1000 (25x40)	Yes	Star Foot	2.80
P60110	1.8	12.3	46.2	1000 (25x40)	No	Round	2.80

TECHNICAL DATA SHEET

Doc No – Ab/TDS/Cryo-27
Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Cryo Vial Internal Threaded Sterile with Silicon Seal, PP

P60112	4.5	12.3	90.3	500 (10x50)	Yes	Star Foot	6.20
P60113	4.5	12.3	89	500 (10x50)	No	Round	6.30

1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.5 Color:

Translucent

1.6 Sterilization:

Sterile catalog number from this lot have been pre-sterilizes and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation. SAL 10^{-6}

1.7 Quality control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/Cryo-27
Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Cryo Vial Internal Threaded Sterile with Silicon Seal, PP

1.8 Other Information:

Recommended for single use only

2. Basic Features

2.1 Temperature range: Tube : -196 °C to + 121 °C

* Caution: Do not use Abdos Cryogenic Vial in the liquid phase of liquid nitrogen. Improper use may cause liquified nitrogen to be trapped inside the vial and lead to pressure build up resulting in possible explosion or biohazard release. Can be used in vapour phase of liquid nitrogen.

2.2 Autoclavability: Autoclavable at 15 psi and 121 °C for 20 minutes (Caps & Vials)

2.3 Leak Test at 635 mm Hg with Cap fitted: No leakage

2.4 Large White frosted portion for easy writing & silk screen accurate graduation for sample identification.

2.5 These Cryo Vials have silicon washer 'O' ring for secure seal.

2.6 Star Foot Bottom base helps in single hand operation.

2.7 Secure tamper evident packaging along with self-stand Zipper pouch.

2.8 No decolorization after sterilizing.

2.9 Manufactured in Clean room ISO Class 8 (100,000)

3.0 Chemical resistance: see website (www.abdoslifesciences.com)

3.1 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/Cryo-27
Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Cryo Vial Internal Threaded Sterile with Silicon Seal, PP

Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

Instruction for Use

Product can be used for food contact also.
Storage temp. Max 30°C
Before using the product check crackness, pouch punctured etc.
Check all the information printed on pouches or boxes.
Always wear gloves while handling specimens.


WARNING

DO NOT store strong oxidizing agents.
DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Consult local fire codes prior to storage of flammable liquids.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/Cryo-27
Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Cryo Vial Internal Threaded Sterile with Silicon Seal, PP

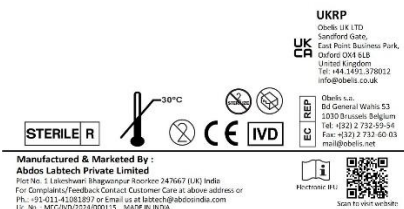

Quality Manager
Abdos Labtech Private Limited
Employee Code 30110

