

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-101

Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.

AmpQlear™ PCR Strip Tubes without Cap 8 Strips, PP

Valid for Cat. No. P10706, P10708, P10732, P10733

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

Strip Tubes without Cap 8 Strips, PP, DNase, RNase, Human DNA, ATP, Pyrogen Free, & PCR Inhibitors Free, IVD, CE & ISO18385 certified.

1.2 Design: Strip Tubes without Cap & with Cap

1.3 Performance Specification:

Code	Capacity (ml)	Type	Height (mm)	No. per Case
P10706	0.2	8 Strip Standard Tubes	20.70	125 Strips
P10708	0.1	8 Strip Low profile Tubes	16.10	125 Strips
P10732	0.2	White 8 Strip Standard Tubes	20.70	125 Strips
P10733	0.1	White 8 Strip Low Profile Tubes	16.10	125 Strips

1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-101

Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.

AmpQlear™ PCR Strip Tubes without Cap 8 Strips, PP

1.5 Color: Translucent (Natural)/White

1.6 Quality control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release

1.7 Other information

1.7.1 DNase, RNase, DNA, Endotoxin and PCR inhibitor free.

1.7.2 Recommended for Single use only

1.7.3 Thin and uniform wall thickness to ensure optimal heat transfer reducing insulating effect.

1.7.4 Snap cap ensures tight fit eliminating risk of sample evaporation. Flat caps are optically clear for real time qPCR.

1.8 Basic Features

1.8.1 Temperature range: -80 °C to + 121°C

1.8.2 Autoclavability: Yes, Autoclavable at 15 psi and 121 °C for 20 minutes

1.8.3 Manufactured in Aseptic Environment

1.8.4 Chemical resistance: see website (www.abdoslifesciences.com)

1.8.5 Tubes have ultra thin walls to reduce the insulting effect on the strips

1.8.6 Strips are available in both Standard & Low Profile wells

1.8.7 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-101

Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.

AmpQlear™ PCR Strip Tubes without Cap 8 Strips, PP

Autoclaving Guideline for Plastics

All plastic consumables are Autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use dry cycle as this may cause deformation or weakening of plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimens.


WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.

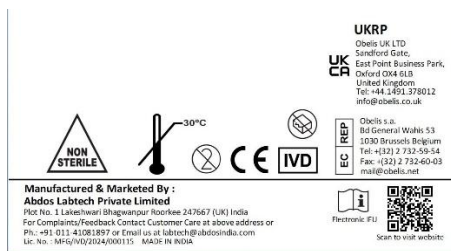

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-101

Rev No: 08 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.

AmpQlear™ PCR Strip Tubes without Cap 8 Strips, PP

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110

