

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-100

Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
AmpQlear™ PCR Tubes, PP

Valid for Cat. No. P10701, P10702

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

1.1.1 PCR Tubes, PP, DNase, RNase, Human DNA, ATP & Pyrogen Free, IVD & ISO 18385 Certified

1.1.2 PCR tubes are made of high purity virgin USP Class VI Medical grade polypropylene (PP) conforming to US FDA 21 CFR, free from natural rubber and heavy metal.

1.1.3 Thin and uniform wall thickness to ensure optimal heat transfer reducing insulating effect.

1.1.4 Attached snap cap with lid. Caps available are in domed & flat shape. Snap cap ensures tight fit eliminating risk of sample evaporation. Flat caps are optically clear for real time qPCR.

1.2 Design:

1.2.1 Tube: Cylindrical upper and conical bottom

1.2.2 Tube cap: Attached Snap Cap with lid

1.2.3 Caps are optically clear for real time qPCR

1.3 Performance Specification:

Code	Capacity (ml)	Type	Height (mm)	No. per Case
P10701	0.2	Clear Flat Cap	21.50	1000(2x500)
P10702	0.2	Domed Cap	23.16	1000(2x500)

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-100

Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
AmpQlear™ PCR Tubes, PP

1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.5 Color: Translucent (Natural)

1.6 Quality control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release

1.7 Other information

1.7.1 DNase, RNase, DNA, Endotoxin and PCR inhibitor free.

1.7.2 Recommended for Single use only

1.7.3 Snap cap ensures tight fit eliminating risk of sample evaporation. Flat caps are optically clear for real time qPCR.

1.8 Basic Features

1.8.1 Temperature range: -80 °C to + 121°C

1.8.2 Autoclavability: Yes, Autoclavable at 15 psi and 121 °C for 20 minutes

1.8.3 Manufactured in Aseptic Environment

1.8.4 Chemical resistance: see website (www.abdoslifesciences.com)

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-100

Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
AmpQlear™ PCR Tubes, PP

1.8.5 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

Autoclaving Guideline for Plastics

All plastic consumables are Autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use dry cycle as this may cause deformation or weakening of plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimens.

WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure


Consult local fire codes prior to storage of flammable liquids.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ PCR-100

Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
AmpQlear™ PCR Tubes, PP

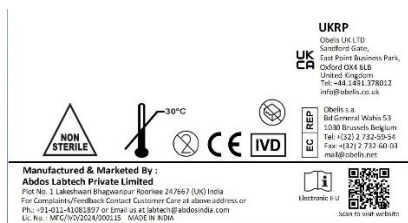

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



t.

