

# Technical Data Sheet

Doc No – Ab/TDS/Cryo-27

Rev No: 05 Date 04/2023-24

## TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Cryo Vial External Threaded Sterile, PP

Valid for Cat. No. P60115, P60116, P60117, P60118, P60116A, P60122, P60123

SRN (Single Registration Number): IN-MF-000022272

### 1. Description/Specification

**1.1 Description:** Disposable Cryo vials, DNase, Rnase, Human DNA, ATP & Pyrogen Free, IVD certified

#### 1.2 Design:

1.2.1 Vials : Available in Star Foot, Round Bottom with External Thread

1.2.3 Cap: Assembled

#### 1.3 Performance Specifications:

1.3.1 Vial Volume: 1.0 ml, 1.8 ml and 4.5 ( ± 0.1 ml)

##### 1.3.1 Graduation and writing surface:

1.3.1.1 Large Printed Graduation

1.3.1.2 White Patch for easy writing on the tube

Code	Capacity (ml)	Outer Diameter (mm)	Height	No. Per Case	Self-Standing	Bottom	Brimful Capacity (Approx ml)
P60115	1.0	12.00	41.00	1000 (25x40)	Yes	Star Foot	1.70
P60116	1.8	12.00	48.5	1000 (25x40)	Yes	Star Foot	2.50
P60117	1.8	12.00	47.5	1000 (25x40)	No	Round	2.50
P60118	4.5	11.80	86.5	500 (10x50)	Yes	Star Foot	5.30
P60116A	1.8	12.00	48.5	1000 (25x40)	Yes	Star Foot	2.50
P60122	4.5	11.80	86.5	500 (10x50)	Yes	Star Foot	5.30
P60123	4.5	11.80	86.5	500 (10x50)	Yes	Round Bottom	5.30

\* Please note compatible coders for P60118 Cryo Vial 4.5 ml Skirted Foot is P60104.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email: labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

## 1.4 Material:

1.4.1 **Tube:** Polypropylene (PP), Hydrophobic **Cap**-HDPE, confirm US FDA 21CFR

1.4.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for Pharmaceutical use.

1.4.3 The product has passed the USP testing including class VI tests, Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE, 2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water

1.4.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore, BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

## 1.5 Color:

Translucent (Natural)

## 1.6 Sterilization:

Sterile catalog number from this lot have been pre-sterilizes and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare Products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation. SAL  $10^{-6}$

## 1.7 Quality Control:

Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

## 1.8 Other Information:

Recommended for single use only

## 1.9. Basic Features

**1.9.1 Temperature range:** Tube:  $-196^{\circ}\text{C}$  to  $+121^{\circ}\text{C}$

\* Caution: Do not use Abdos Cryogenic Vial in the liquid phase of liquid nitrogen. Improper use may cause liquified nitrogen to be trapped inside the vial and lead to pressure build up resulting in possible explosion or biohazard release. Can be used in vapour phase of liquid nitrogen.

**1.9.2 Autoclavability:** Autoclavable at 15 psi and  $121^{\circ}\text{C}$  for 20 minutes (Caps & Vials)

**1.9.3 Leak Test at 635 mm Hg with Cap fitted:** No leakage

**1.9.4 Large White frosted portion** for easy writing & silk screen accurate graduation for sample identification.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email: labtech@abdosindia.com  
Web: [www.abdoslivesciences.com](http://www.abdoslivesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

**1.9.5** These cryo vial are designed for mechanical freezing and screw cap reduces risk of contamination due to handling.

**1.9.6** Star Foot Bottom base helps in single hand operation.

**1.9.7** Secure tamper evident packaging along with self-stand Zipper pouch.

**1.9.8** No discolorization after re-sterilizing.

**1.9.9** Manufactured in Clean room ISO Class 8 (100,000)

**2.0 Chemical resistance:** see website ([www.abdoslifesciences.com](http://www.abdoslifesciences.com))

## **Autoclaving Guideline for Plastics**

All plastic consumables are autoclavable at 121<sup>0</sup>C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50<sup>0</sup>C before removing product.

## **Instruction for Use**

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.

## **WARNING**

**DO NOT** store strong oxidizing agents.

**DO NOT** place any plastic labware in a flame.

**DO NOT** mix any chemicals that may result in a thermic reaction, which can cause product failure


Consult local fire codes prior to storage of flammable liquids.

# Technical Data Sheet

UKRP

Obelis UK LTD  
Sandford Gate,  
East Point Business Park,  
Oxford OX4 6LB  
United Kingdom  
Tel: +44.1491.378012  
info@obelis.co.uk

STERILE R

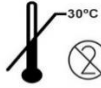

30°C



EC REP



Obelis s.a.  
Bd General Wahis 53  
1030 Brussels Belgium  
Tel: +(32) 2 732-59-54  
Fax: +(32) 2 732-60-03  
mail@obelis.net



**Abdos Labtech Pvt. Ltd.**

Plot No. 1 Lakeshwari Bhagwanpur Roorkee 247667 (UK) India  
Email : labtech@abdosindia.com, www.abdoslifesciences.com  
Ph: +(91) 11 41081897



Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



**Abdos Labtech Private Limited**  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakhewari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903

**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

**Ph. No: +91- (0) 1332235277**  
**Email:labtech@abdosindia.com**  
**Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)**  
**GST: 05AAGCA8605A1Z7**