

Technical Data Sheet

Doc No– Ab/TDS/Micro-83-V

Rev No: 07 Date 05/2024-25

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Petri Dish, (PS) 3 Compartment Aseptic/Sterile

Valid for Cat. No. P10919, P10920, P10921, P10922, P10923, P10924, P10925, P10926

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Petri Dish 3 compartment- Gama Radiated Sterile/Aseptic, DNase, RNase, Human DNA, ATP, Pyrogen Free, CE & IVD certified, ISO24998

1.2 Design: Shallow dish design

1.3 Performance Specifications:

Code	Sterile	Diameter (mm)	Height (mm)	Compartment	Vented	Packing	No. Per Case
P10919	Aseptic	90	15.80	3	3 Vent	Sleeve of 20	480
P10920	Aseptic	90	15.80	3	3 Vent	Individual Wrapped	380
P10921	Aseptic	90	15.80	3	0 Vent	Sleeve of 20	480
P10922	Aseptic	90	15.80	3	0 Vent	Individual Wrapped	380
P10923	Sterile	90	15.80	3	3 Vent	Sleeve of 20	480
P10924	Sterile	90	15.80	3	3 Vent	Individual Wrapped	380
P10925	Sterile	90	15.80	3	0 Vent	Sleeve of 20	480
P10926	Sterile	90	15.80	3	0 Vent	Individual Wrapped	380

1.4 Material:

1.4.1 Petri Dish made up of Hydrophobic Polystyrene (PS), Confirm with US FDA 21CFR

1.4.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for pharmaceutical use.

Technical Data Sheet

1.4.3 The product has passed the USP testing including class VI tests, Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE, 2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe to use in contact with food stuff, Pharmaceuticals & drinking water

1.4.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.5 Color: Transparent (Natural color)

1.6 Quality control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.7 Other information

1.7.1 Petri Dishes are suitable for all bacteriological applications.

1.7.2 Vented dishes have an elevated lid that allows for superior gas exchange and evaporation.

1.7.3 Petri dishes retain their shape up to 55°C.

1.7.4 Easy to Stack, 12° inclination 10 stacking filed with 20 ml media

1.7.5 Recommended for single use only

1.8 Basic Features

1.8.1 Temperature range: Stable up to 60°C.

1.8.2 Autoclavability: NO

1.8.3 Chemical resistance: See website (www.abdoslifesciences.com)

1.8.4 Manufactured in Clean room ISO Class 8 (100,000)

1.8.5 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

1.9 Sterilization:

Sterile catalog number from this lot have been pre-sterilized and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation. SAL 10^{-6}


Technical Data Sheet

Instruction for Use

Product can be used for food contact also.
Storage temp. Max 30°C
Before using the product check crackness, pouch punctured etc.
Check all the information printed on pouches or boxes.
Always wear gloves while handling specimens.

WARNING

DO NOT store strong oxidizing agents.
DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Consult local fire codes prior to storage of flammable liquids.

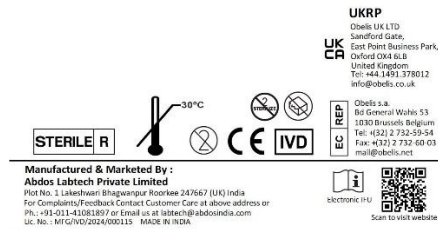

Quality Manager
Abdos Labtech Private Limited
Employee Code 30110

ABDOS LABTECH PVT. LTD.
QUALITY
DEPT.
ROORKEE, INDIA

