

Technical Data Sheet

Doc No – Ab/TDS/LH-65

Rev No : 02 Date 02/2022-23

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

PureRight Last Drop Filter Pipette Tips, Individually Wrapped, PP

Valid for Cat. No. P11040, P11041, P11042, P11043, P11044, P11045, P11046, P11070, P11071

1. Description/Specification

1.1 Description: Disposable PP Micropipette Tips, DNase, RNase, Human DNA, ATP & Pyrogen Free, CE & IVD certified.

1.2 Design:

1.2.1 Open Cylindrical Upper and conical Apex with orifice for liquid dispensing. Precise graduation mark for error free pipetting.

1.3 Performance Specifications:

| Code | Capacity | Colour | No. per Case |
|---------------|-------------------------------|---------|--------------|
| P11040 | 0.1-10µl (Graduated) | Natural | 350 |
| P11041 | 0.1-10µl Extended (Graduated) | Natural | 350 |
| P11042 | 1-20µl (Graduated) | Natural | 350 |
| P11043 | 1-100µl (Graduated) | Natural | 350 |
| P11044 | 1-200µl (Graduated) | Natural | 350 |
| P11045 | 1 – 300µl (Graduated) | Natural | 350 |
| P11046 | 100-1000ul | Natural | 350 |
| P11070 | 100-1000ul XL | Natural | 100 |
| P11071 | 100-1250ul | Natural | 100 |

1.3.2 Fitting to pipettes: Universal Fit (For detail refer the Abdos compatibility chart)

1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 The product has passed the USP testing including class VI tests ,Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE ,2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email: labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.4.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.5 Color:

Natural

1.6 Sterilization:

Sterile catalog number from this lot have been pre-sterilizes and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation. SAL 10^{-6}

1.7 Quality control:

Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.8 Other information:

The retention level is minimal, diamond polished moulds for an extremely hydrophobic & smooth surface to ensure lowest retention in its class. Super filter provides aerosol & aqueous liquid by pass prevention in over pipetting thereby preventing cross contamination. Hydrophobic filter that is additive free & does not absorb samples.

Tips are of single used only.

1.9 Basic Features:

1.9.1. Made up of : High Purity Virgin Polypropylene (PP)

1.9.2. Temperature range : -80°C to $+121^{\circ}\text{C}$

1.9.3. Autoclavability: Can be autoclave but might loose the fitting of filters.

1.9.4. Chemical resistance: see website (www.abdoslifesciences.com)

1.9.5. Retention: Best in its class. (With Acidic & Basic dyes)

1.9.6 Manufactured in Clean room ISO Class 8 (100,000)

1.9.7 ABDOS Super Filter Tips has a **100% pure, porous polyethylene (PE) filter** that prevents aerosol contamination and carryover generated during sample aspiration. In addition, Filter Pipette Tip does not trap samples and allows for sample recovery in the event of over-pipetting. 99.9% Bacterial Aerosol Filtration Efficiency (BFE), (ASTM F21012) Inert, hydrophobic PE filter does not trap or block the liquid. Free of leachables, extractables and heavy metals. Non Cytotoxic (ISO 10993-5 MRM), Minimal Essential Media Elution, GLP Compliant Non-hemolytic.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email: labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7



Technical Data Sheet

Instruction for Use

Product can be used for food contact also.
Storage temp. Max 30°C
Before using the product check crackness, pouch punctured etc.
Check all the information printed on pouches or boxes.
Always wear gloves while handling specimen.

WARNING

DO NOT store strong oxidizing agents.
DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Consult local fire codes prior to storage of flammable liquids.

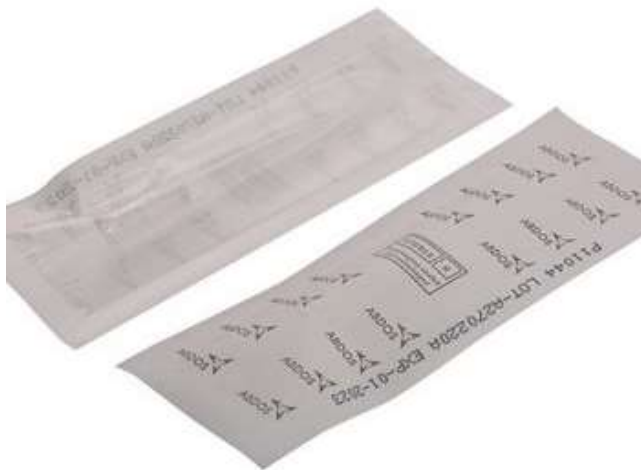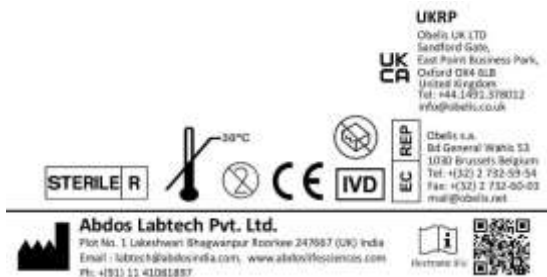

Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakshwari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7