

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Narrow Mouth Bottle, LDPE

**Valid for Catalog No. P11141, P11142, P11143, P11144, P11145, P11146, P11147,
P11148, P11149**

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

Narrow Mouth Bottle, LDPE

- 1.1.1 These bottles come with inner thread for superior sealing and better security.
- 1.1.2 These bottles are used for packing, shipping and storage purposes.
- 1.1.3 Superior chemical resistance.
- 1.1.4 Complete emptying is guaranteed through inclination of shoulders.
- 1.1.5 Suitable for storage of low-density liquids.

1.2 Raw material:

1.2.1 Body: LDPE (Low Density polyethylene)

1.2.2 Cap: PP (Polypropylene)

1.3 Material:

- 1.3.1 Low density polyethylene (LDPE), Hydrophobic, confirm US FDA 21CFR
- 1.3.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.
- 1.3.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.4 Quality Control:

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Narrow Mouth Bottle, LDPE

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (ml)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11141	60	17.5	82.4	72
P11142	125	21	97.6	72
P11143	250	21	130.7	72
P11144	500	24.6	165.5	48
P11145	4	12.2	39	72
P11146	8	17.5	42.3	72
P11147	15	17.5	56.1	72
P11148	30	17.5	60	72
P11149	1000	52.1	221.5	24

1.6 GENERAL FEATURES:

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Narrow Mouth Bottle, LDPE

CE	Yes
Autoclavability	No
Translucence	Yes
Resistant	To be affected by radiations
Sterilization	NA
Temperature	From – 40°C to + 80°C
Leak Proof	Yes

1.7 Chemical resistance: See website (www.abdoslifesciences.com)

1.7.1 Manufactured in Aseptic Environment

1.7.2 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.

1.8.3 Poor resistance; exposure to chemicals is not allowed, or causes immediate permanent damage.

Results of testing with reference materials are to be considered as indication: in case of specific use, it is recommended to carry out preliminary testing.

Information relating to the resistance of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Narrow Mouth Bottle, LDPE

WARNING

DO NOT store strong oxidizing agents.
DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.
Consult OSHA guidelines prior to handling or lifting filled carboys/bottles/aspirators.

If applicable: Using the spigot as a handle can ruin the threads and cause the spigot not to seal properly.

Cleaning:

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

Cap and Spigot Cleaning: (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.





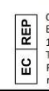

Check all the information printed on pouches or boxes.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Narrow Mouth Bottle, LDPE

Always wear gloves while handling specimens.

 NON STERILE	 30°C	 CE	 IVD	 UK REP	UKRP Obelis UK LTD Sandford Gate, East Point Business Park, Oxford OX4 6LB United Kingdom Tel: +44.1491.378012 info@obelis.co.uk
Manufactured & Marketed By : Abdos Labtech Private Limited Plot No. 1, Lakshwari Bhagwanpur Roorkee 247667 (UK) India For Complaints/Feedback Contact Customer Care at above address or Ph.: +91-011-41081897 or Email us at labtech@abdosindia.com Lic. No.: MFG/IVD/2024/000115 MADE IN INDIA					 Electronic IFU Scan to visit website

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-138
Rev No: 07 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Narrow Mouth Bottle, LDPE

