

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Carboy with Stopcock, PP

Valid for Catalog No. P11520, P11521, P11524

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

Carboy with Stopcock, PP

1.1.1 Abdos Carboys are manufactured using high purity virgin USP Class VI Medical grade polypropylene (PP).

1.1.2 These carboys come with secure sealing cap for better security and long-term storage.

1.1.3 These carboys are used for storing and shipping reagents, chemicals, culture media, and standards.

1.1.4 These are very strong and durable for extended use over time with a strong carrying handle on both sides of the carboy.

1.1.5 Comes with high quality and durable stopcock or tap for easy withdrawal of liquids, chemicals, reagents, culture media and standards

1.1.5 Can be used at a temperature below 0°C.

1.1.6 Superior chemical resistance.

1.2 Raw material:

1.2.1 Body: PP (Poly Propylene)

1.2.2 Cap: PP (Polypropylene)

1.3 Material:

1.3.1 Poly Propylene (PP) & High Density poly Ethylene (HDPE), Hydrophobic, confirm US FDA 21CFR

1.3.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Carboy with Stopcock, PP

1.3.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.4 Quality Control:

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (Ltr)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11520	10	63.8	378	1
P11521	20	64.2	525	1
P11524	Spare Stopcock			1

1.6 GENERAL FEATURES:

CE	Yes
Autoclavability	Yes
Translucence	Yes
Resistant	To be effected by radiations

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Carboy with Stopcock, PP

Sterilization	NA
Temperature	From 0°C to 121°C
Leak Proof	Yes

1.7 Chemical resistance: See website (www.abdoslivesciences.com)

1.7.1 Manufactured in Aseptic Environment

1.7.2 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.

1.8.3 Poor resistance; exposure to chemicals is not allowed or causes immediate permanent damage.

Results of testing with reference materials are to be considered as indication: in case of specific use, it is recommended to carry out preliminary testing.

Information relating to the resistance of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products.

WARNING

DO NOT store strong oxidizing agents.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Carboy with Stopcock, PP

DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.
Consult OSHA guidelines prior to handling or lifting filled carboys/bottles/aspirators.

If applicable: Using the spigot as a handle can ruin the threads and cause the spigot not to seal properly.

Cleaning:

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

Cap and Spigot Cleaning: (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.


Check all the information printed on pouches or boxes.

Always wear gloves while handling specimens.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Carboy with Stopcock, PP

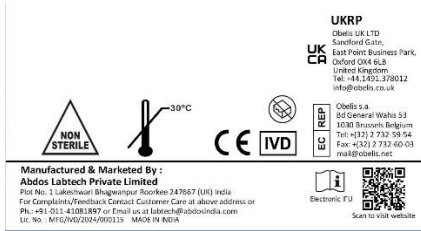

Quality Manager
Abdos Labtech Private Limited
Employee Code 30110



TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-148
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Carboy with Stopcock, PP

