

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-152
Rev No: 06 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Dropping Bottle Euro Type, LDPE

Valid for Catalog No. P11305, P11306, P11307, P11308

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Dropping Bottle Euro Type, LDPE

1.1.1 These dropping bottles are used for routine storage and pouring applications in laboratory.

1.1.2 These bottles can be squeezed to ensure precise and steady dropping of liquid from the bottle.

1.1.3 Moulded from high quality virgin low-density polyethylene to ensure excellent results.

1.1.4 Superior chemical resistance.

1.2 Raw material:

1.2.1 Body: LDPE (Low Density poly Ethylene)

1.2.2 Cap: PP (Polypropylene)

1.3 Material:

1.3.1 Low density polyethylene (LDPE), Hydrophobic, confirm US FDA 21CFR

1.3.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.3.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.4 Quality Control:

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-152
Rev No: 06 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Dropping Bottle Euro Type, LDPE

Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (ml)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11305	30	17.5	60	12
P11306	60	17.5	82.4	12
P11307	125	23.6	94.2	12
P11308	205	22.8	101.8	12

1.6 GENERAL FEATURES:

CE	Yes
Autoclavability	No
Translucence	Yes
Resistant	To be effected by radiations
Sterilization	NA
Temperature	From – 40°C to + 80°C
Leak Proof	Yes

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-152
Rev No: 06 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Dropping Bottle Euro Type, LDPE

1.7 Chemical resistance: See website (www.abdoslifesciences.com)

1.7.1 Manufactured in Aseptic Environment

1.7.2 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.

1.8.3 Poor resistance; exposure to chemical is not allowed or causes immediate permanent damage.

Results of testing with reference materials are to be considered as indication: in case of specific use, it is recommended to carry out preliminary testing.

Information relating to the resistances of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products.

WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-152
Rev No: 06 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Dropping Bottle Euro Type, LDPE

If applicable: Using the spigot as a handle can ruin the threads and causing the spigot not to seal properly.

Cleaning:

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

Cap and Spigot Cleaning: (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

Instruction for Use


Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.

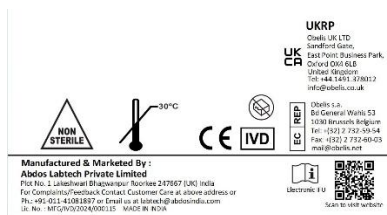

TECHNICAL DATA SHEET

Doc No – Ab/TDS/BC-152
Rev No: 06 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Dropping Bottle Euro Type, LDPE

Quality Manager
Abdos Labtech Private Limited
Employee Code 30110

