

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Racked Centrifuge Tube Conical Bottom, Paper

Valid for Cat. No. Paper Rack: P10409 (15ml), P10410 (50ml)

Racked Non-Sterile: P10411, P10413, P10415, P10417

Racked Sterile: P10412, P10414, P10416, P10418, P10462, P10463, P10455

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Disposable PP/PS Tube, HDPE Cap & Paper rack, DNase, RNase, Human DNA, ATP & Pyrogen Free, CE & IVD certified

1.2 Design:

1.2.1 Tube: Cylindrical Upper and conical Bottom

1.2.2 Tube Cap: Separate Threaded HDPE Cap

1.2.3 Rack: These racks are made of eco-friendly paper to hold 25 Abdos 15 ml & 50 ml Conical Bottom Centrifuge Tubes.

1.3 Performance Specifications:

Code	Capacity (ml)	OD (mm)	Height (mm)	Sterile	Tube	Cap Colour	No./Case
P10409	15	NA	47	NO	NA	White	10
P10410	50	NA	48	NO	NA	White	10
P10411	15	16.5	121	NO	Natural	Blue	500(25pcsx20racks)
P10412	15	16.5	121	YES	Natural	Blue	500(25pcsx20racks)
P10462	15	16.5	121	YES	Natural	White	500(25pcsx20racks)

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Racked Centrifuge Tube Conical Bottom, Paper

P10413	15	16.5	121	NO	Amber	Blue	500(25pcsx20racks)
P10414	15	16.5	121	YES	Amber	Blue	500(25pcsx20racks)
P10455	15	16.5	121	YES	Natural	Red	500(25pcsx20racks)
P10415	50	29.2	116	NO	Natural	Blue	500(25pcsx20racks)
P10416	50	29.2	116	YES	Natural	Blue	500(25pcsx20racks)
P10417	50	29.2	116	NO	Amber	Blue	500(25pcsx20racks)
P10418	50	29.2	116	YES	Amber	Blue	500(25pcsx20racks)
P10463	50	29.2	116	YES	Natural	White	500(25pcsx20racks)

1.4 Material:

Rack- Paper

1.4.1 Tube - Polypropylene (PP) & Cap- High Density Polyethylene (HDPE) confirm US FDA 21CFR

1.4.2 The product complies with USP class VI. Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Racked Centrifuge Tube Conical Bottom, Paper

1.5 Color:

Tube color- Translucent (Natural), and Amber

Cap color- Blue/White/Red

Paper Racked Color: White

1.6 Graduation and writing surface:

1.6.1 Large Printed Graduation

1.6.2 White Patch for easy writing on the tube

1.6.3 Printed Max Freeze Level Indicator with Red Colour

1.6.5 Reverse Graduation

1.7 Sterilization:

Sterile catalog number from this lot have been pre-sterilizes and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation ,SAL 10^{-6}

1.8 Quality Control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.9 Other Information:

Recommended for single use only

Pack Size: P10409 & P10410- 10 pcs/case

P10411, P10413, P10415, P10417, P10463– **500 Pcs (25X 20 Racks) case**

P10412, P10414, P10416, P10418, P10462, P10455– **500 Pcs (25X 20 Racks) case**

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Racked Centrifuge Tube Conical Bottom, Paper

2. Basic Features

- 2.1 Temperature range:** Tube: -80⁰C to +121⁰C (PP), -10⁰C to + 70⁰C (PS)
- 2.2 Autoclavability:** Only PP is Autoclavable at 15 psi and 121⁰C for 20 minutes
- 2.3 Centrifugation:** 18,000 RCF (Rated with Double Distilled water at 25⁰C for 30 minutes)
(Only PP Tubes)
- 2.4 Leak Test at 635 mm Hg:** No leakage
- 2.5** These racks are made of eco-friendly paper to hold 25 pieces of CT 15 & 50ml.
- 2.6** Manufactured in Clean room ISO Class 8 (100,000)
- 2.7 Chemical resistance:** see website (www.abdoslifesciences.com)
- 2.8 Intended use:** Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Racked Centrifuge Tube Conical Bottom, Paper

Instruction for Use

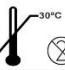




Product can be used for food contact also.
Storage temp. Max 30°C
Before using the product check crackness, pouch punctured etc.
Check all the information printed on pouches or boxes.
Always wear gloves while handling specimen.
Ensure all sample tubes are evenly filled
For each tube inserted in the rotor, add a tube of equal weight directly opposite it.

WARNING

DO NOT store strong oxidizing agents.
DO NOT place any plastic labware in a flame.
DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure
Wear a face shield and / or safety goggles if you have to work anywhere near a centrifuge that's in use.
Do not bump, jar, or move the centrifuge while the rotor is spinning. Make sure you don't have the cord dangling from a table edge where someone could catch their foot in it and pull down the centrifuge
Consult local fire codes prior to storage of flammable liquids.


UKRP
Obelis UK LTD
Sandford Gate,
East Point Business Park,
Oxford OX4 6LB
United Kingdom
Tel: +44 1851 378012
info@obelis.co.uk

UK
CE
CR

STERILE R     

Obelis S.A.
Bd General Wafes 53
2030 Brussels Belgium
Tel: +32 2 732 55 54
Fax: +32 2 732 60 03
mail@obelis.net

Manufactured & Marketed By:
Abdos Labtech Private Limited
Plot No. 1, Lakhewari Bhagwanpur Roorkee 247667 (UK) India
For Complaints/Feedback Contact Customer Care at above address or
Ph: +91-011-41083897 or Email us at labtech@abdosindia.com
Lic. No.: MFG/IVD/2024/000155 MADE IN INDIA

 **electronic IVD**
Scan to visit website

TECHNICAL DATA SHEET

Doc No – Ab/TDS/CT-09
Rev No: 10 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Racked Centrifuge Tube Conical Bottom, Paper

Quality Manager
Abdos Labtech Private Limited
Employee Code 30110

