

# TECHNICAL DATA SHEET

Doc No – Ab/TDS/ Micro-85

Rev No: 08 Date- 05/2024-25

**Abdos Labtech Pvt. Ltd.**

**Soft Inoculation Loop Gama Radiated Sterile, HIPS**

Valid for Cat. No. P10985, P10986, P10987, P10988, P10977, P10978, P10979, P10980

**SRN (Single Registration Number):** IN-MF-000022272

## 1. Description/Specification

### 1.1 Description: Soft Inoculation Loop Gama Radiated Sterile, HIPS

1.1.1 Smooth surface provides uniform and gentle streak without damaging the gel surface.

1.1.2 Round shaft improves grip, assists orientation and makes operation easy.

1.1.3 Available in individual and zip lock sleeve that eliminates risk of contamination.

### 1.2 Design: Round shaft Loop

### 1.3 Performance Specifications:

| Code          | Capacity (µl)                   | Colour | Length (mm) | Packing            | No. Per Case |
|---------------|---------------------------------|--------|-------------|--------------------|--------------|
| <b>P10985</b> | 1µl loop + Incorporated Needle  | White  | 200         | Sleeve of 10       | 500          |
| <b>P10986</b> | 1µl loop + Incorporated Needle  | White  | 200         | Individual Wrapped | 500          |
| <b>P10987</b> | 10µl loop + Incorporated Needle | Blue   | 200         | Sleeve of 10       | 500          |
| <b>P10988</b> | 10µl loop + Incorporated Needle | Blue   | 200         | Individual Wrapped | 500          |
| <b>P10977</b> | 1µl loop + Blunt                | White  | 200         | Sleeve of 10       | 500          |

# TECHNICAL DATA SHEET

Doc No – Ab/TDS/ Micro-85

Rev No: 08 Date- 05/2024-25

**Abdos Labtech Pvt. Ltd.**

**Soft Inoculation Loop Gama Radiated Sterile, HIPS**

|               |                   |       |     |                       |     |
|---------------|-------------------|-------|-----|-----------------------|-----|
| <b>P10978</b> | 1µl loop + Blunt  | White | 200 | Individual<br>Wrapped | 500 |
| <b>P10979</b> | 10µl loop + Blunt | Blue  | 200 | Sleeve of 10          | 500 |
| <b>P10980</b> | 10µl loop + Blunt | Blue  | 200 | Individual<br>Wrapped | 500 |

## 1.4 Material:

1.4.1 Inoculation Loop made up of HIPS (High Impact Polystyrene)

1.4.2 Confirm with US FDA Title 21CFR of US Food and Drug Administration

1.4.3 Free from Heavy metals, It is safe use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.4 Product is manufactured without the use of Raw Materials of animal origin Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

**1.5 Color:** White & Blue

## 1.6 Quality control:

Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

# TECHNICAL DATA SHEET

Doc No – Ab/TDS/ Micro-85

Rev No: 08 Date- 05/2024-25

**Abdos Labtech Pvt. Ltd.**  
**Soft Inoculation Loop Gama Radiated Sterile, HIPS**

---

## 1.7 Other information

1.7.1 Impact resistance and tough.

1.7.2 Superior for its hardness, gloss, toughness, and electrical insulation properties.

## 1.8 Basic Features

**1.8.1 Temperature range:** -10°C to 70°C

**1.8.2 Autoclavability:** NO

**1.8.3 Chemical resistance:** See website. ([www.abdoslifesciences.com](http://www.abdoslifesciences.com))

**1.8.4 Manufactured in Clean room ISO Class 8 (100,000)**

**1.8.5 Intended use:** Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

## 1.9 Sterilization:

**Sterile catalog number** from this lot have been pre-sterilizes and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products Requirements for validation and Routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation. SAL  $10^{-6}$

### Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimens.

# TECHNICAL DATA SHEET

Doc No – Ab/TDS/ Micro-85

Rev No: 08 Date- 05/2024-25

**Abdos Labtech Pvt. Ltd.**

**Soft Inoculation Loop Gama Radiated Sterile, HIPS**


## WARNING

**DO NOT** store strong oxidizing agents.

**DO NOT** place any plastic labware in a flame.

**DO NOT** mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.

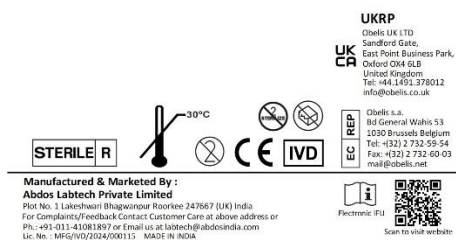

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



# TECHNICAL DATA SHEET

Doc No – Ab/TDS/ Micro-85

Rev No: 08 Date- 05/2024-25

**Abdos Labtech Pvt. Ltd.**

**Soft Inoculation Loop Gama Radiated Sterile, HIPS**

