

Technical Data Sheet

Doc No–Ab/TDS/BC-152

Rev No: 03 Date 04/2023-24

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

Carboy Wide Mouth, PP

Valid for Catalog No. P11522, P11523

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description:

Carboy Wide Mouth, PP

1.1.1 Abdos Carboys are manufactured using high purity virgin USP Class VI Medical grade polypropylene (PP).

1.1.2 These carboys come with secure sealing cap for better security and long-term storage.

1.1.3 These carboys are used for storing and shipping reagents, chemicals, culture media, and standards.

1.1.4 These are very strong and durable for extended use over time with a strong carrying handle on both sides of the carboy.

1.1.5 Comes with wide mouth easy withdrawal and refill of liquids, chemicals, reagents, culture media and standards.

1.1.6 Can be used at a temperature below 0°C.

1.1.7 Superior chemical resistance.

1.2 Raw material:

1.2.1 Body: PP (Poly Propylene)

1.2.2 Cap: HDPE (High Density poly Ethylene)

1.3 Material:

1.3.1 Poly Propylene (PP), Hydrophobic, confirm US FDA 21CFR & High Density poly Ethylene (HDPE), Hydrophobic, confirm US FDA 21CFR

1.3.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for pharmaceutical use.

1.3.3 The product has passed the USP testing including class VI tests, Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.3.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE, 2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & its following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email: labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.3.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

1.4 Quality Control:

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (Ltr)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11522	10	76	378	1
P11523	20	85	525	1

1.6 GENERAL FEATURES:

CE	Yes
Autoclavability	Yes
Translucence	Yes
Resistant	To be effected by radiations
Sterilization	NA
Temperature	From 0°C to 121°C
Leak Proof	Yes

1.7 Chemical resistance: See website (www.abdoslifesciences.com)

1.7.1 Manufactured in Aseptic Environment

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.

1.8.3 Poor resistance; exposure to chemical is not allowed or causes immediate permanent damage. Results of testing with reference materials are to be considered as indication: in case of specific use, it is recommended to carry out preliminary testing.

Information relating to the resistances of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products.

WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.

Consult OSHA guidelines prior to handling or lifting filled carboys/bottles/aspirators.

If applicable: Using the spigot as a handle can ruin the threads and causing the spigot not to seal properly.

Cleaning:

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

Cap and Spigot Cleaning: (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakheswari, Bhagwanpur,
Roorkee-247667, (UK) INDIA


Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903
Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7

Technical Data Sheet

UKRP

Obelis UK LTD
Sandford Gate,
East Point Business Park,
Oxford OX4 6LB
United Kingdom
Tel: +44.1491.378012
info@obelis.co.uk

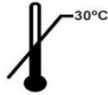

Obelis s.a.
Bd General Wahis 53
1030 Brussels Belgium
Tel: +(32) 2 732-59-54
Fax: +(32) 2 732-60-03
mail@obelis.net



Abdos Labtech Pvt. Ltd.

Plot No. 1 Lakhewari Bhagwanpur Roorkee 247667 (UK) India
Email : labtech@abdosindia.com, www.abdoslifesciences.com
Ph: +(91) 11 41081897



Electronic IFU



Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



Abdos Labtech Private Limited
(A Division of Abdos Group)
(An ISO 9001,13485,14001,18385,GMP & CE certified company)
Shiv Ganga Industrial Estate
Plot No-1, Lakhewari, Bhagwanpur,
Roorkee-247667, (UK) INDIA

Registered office
13A Government
Place East, 2nd Floor Kolkata
West Bengal – 700069
CIN: U85100WB2008PTC123903

Corporate Office
319, 3rd Floor, DLF Tower A, DDA Dist. Centre
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277
Email:labtech@abdosindia.com
Web: www.abdoslifesciences.com
GST: 05AAGCA8605A1Z7