

# Technical Data Sheet

Doc No – Ab/TDS/LH-62

Rev No : 04 Date 02/2022-23

## TECHNICAL DATA SHEET

**Abdos Labtech Pvt. Ltd.**

**Filter Pipette Tips, Racked PP**

**Valid for Cat. No. P11020, P11022, P11024, P11026, P11028, P11032, P11034, P11036, P11038, P11039, P11030, P11031**

### 1. Description/Specification

**1.1 Description:** Disposable PP Micropipette Racked Tips, DNase, RNase, Human DNA, ATP & Pyrogen Free, CE & IVD certified.

#### 1.2 Design:

**Tip:** Open Cylindrical Upper and conical Apex with orifice for liquid dispensing. Precise graduation mark for error free pipetting.

**Box:** Rectangular Box with Transparent Top Cover with slip lock.

Removable Micro tip Tray Suitable to Work in a Robotic System.

Compatible for use with Multichannel Pipettes

Top cover mark A-H in Y-axis for easy discrimination.

#### 1.3 Performance Specifications:

Code	Capacity (µl)	No./Case	Colour
P11020	0.1-10 (Graduated)	10Racksx96pcs	Natural
P11022	0.1-10 Extended (Graduated)	10Racksx96pcs	Natural
P11024	1-20 ( Graduated)	10Racksx96pcs	Natural
P11026	1-50 (Graduated)	10Racksx96pcs	Natural
P11028	1-100 (Graduated)	10Racksx96pcs	Natural
P11032	1-200 (Graduated)	10Racksx96pcs	Natural
P11034	1-300 (Graduated)	10Racksx96pcs	Natural
P11036	100-1000ul	10Racksx96pcs	Natural
P11038	100-1000ul XL	10Racksx96pcs	Natural
P11039	100-1250ul	10Racksx96pcs	Natural
P11030	100-5000ul	6Racksx24pcs	Natural
P11031	500-10000ul	6Racksx24pcs	Natural



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslivesciences.com](http://www.abdoslivesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

1.3.2 Fitting to pipettes: Universal Fit (For detail refer the Abdos compatibility chart)

## 1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 The product has passed the USP testing including class VI tests, Materials has successfully passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE, 2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & its following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water

1.4.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

## 1.5 Color:

**Tip Colour** -Natural

**Wafer Colour** - 10ul Green, 20 to 100ul Yellow, 200 to 300ul Orange, 1000ul Blue, 1000XL & 1250ul Light Blue

## 1.6 Sterilization:

NA

## 1.7 Quality control:

Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

## 1.8 Other information:

Reusable, only box can be autoclavable up to 10 autoclave cycle. Super filter provides aerosol & aqueous liquid by pass prevention in over pipetting thereby preventing cross contamination. Hydrophobic filter that is additive free & does not absorb samples. Gen 2 Racked Filter Tips are offered in prefilled racks that do not require manual accumulation in empty boxes. Polypropylene racks are very robust, sturdy with consistent and thick wall thickness. Gen 2 Racks are hinged with transparent lid & base providing excellent view inside the box. Gen 2 Racks are stackable and reusable once empty.

Tips are of single used only.

## 1.9. Basic Features:

1.9.1 Made up of : High Purity Virgin Polypropylene (PP)

1.9.2 Temperature range : -80°C to + 121°C for 20 minutes



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001, 13485, 14001, 18385, GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email: labtech@abdosindia.com  
Web: [www.abdoslivesciences.com](http://www.abdoslivesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

**1.9.3 Autoclavability:** Can be autoclave but might loose the fitting of filters.

**1.9.4 Chemical resistance:** see website (www.abdoslifesciences.com)

**1.9.5 Retention:** Best in its class. (with Acidic & Basic dyes)

**1.9.6 Manufactured in Clean room ISO Class 8 (100,000)**

**1.9.7 ABDOS Super Filter Tips** has a **100% pure, porous polyethylene (PE) filter** that prevents aerosol contamination and carryover generated during sample aspiration. In addition, Filter Pipette Tip does not trap samples and allows for sample recovery in the event of over-pipetting. 99.9% Bacterial Aerosol Filtration Efficiency (BFE), (ASTM F21012) Inert, hydrophobic PE filter does not trap or block the liquid. Free of leachables, extractables and heavy metals. Non Cytotoxic (ISO 10993-5 MRM , Minimal Essential Media Elution, GLP Compliant Non-hemolytic.

## Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121<sup>0</sup>C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50<sup>0</sup>C before removing product.

## Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.

## WARNING

**DO NOT** store strong oxidizing agents.

**DO NOT** place any plastic labware in a flame.

**DO NOT** mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.



**Abdos Labtech Private Limited**  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakshewari, Bhagwanpur,  
Roorkh-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

**Ph. No: +91- (0) 1332235277**  
**Email:labtech@abdosindia.com**  
**Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)**  
**GST: 05AAGCA8605A1Z7**


# Technical Data Sheet



## Abdos Life Science successfully filed **PATENT** in **Pipette Tip Box Category**

*Anti-Warpage Technology* offers strength & security at liquid handling

Transparent top cover & bottom provides an excellent view of the Micro Tips inside the box

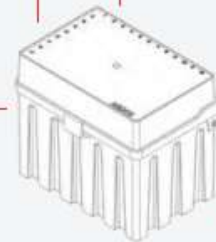

Design with consistent wall thickness for repetitive use in single and multi-channel pipettes

Hinged lid provides fast, safe pipette tip access

Green Manufacturing & Reusable once empty, saves plastics



**Anti-Warpage  
Technology**



**Anti-Shrink  
Technology**



**Robust Design  
& Strength**



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email: labtech@abdosindia.com  
Web: [www.abdoslivesciences.com](http://www.abdoslivesciences.com)  
GST: 05AAGCA8605A1Z7