

Technical Data Sheet

Doc No Ab/TDS/Eqpt-255

Rev No: 02 Date 01/2021

TECHNICAL DATA SHEET

Abdos Labtech Pvt. Ltd.

WAVEMAX – Incu- Shaker Mini

Valid for Catalog No. E11920, E11925

1. Description Specifications:

1.1 Description: WAVEMAX – Incu- Shaker Mini

1.1.1 Wavemax-Incu-Shaker Mini is a compact shaking incubator that can accept a variety of flasks, tubes, trays, etc.

1.1.2 The standard platform features a non-slip surface ideal for staining trays, tissue culture flasks, and Petri dishes.

1.1.3 Abdos Incu-Shaker offer exceptional temperature uniformity and accuracy due to a proprietary heat distribution technology that ensures temperatures are precisely maintained throughout the chamber within 0.25%, even in the corners.

1.1.4 Featuring the MAGIC Clamp platform, this shaker lets you customize the size and types of vessels being used by simply exchanging clamp sizes. Clamps utilize a powerful magnetic base to securely attach to the shaking platform without screws, bolts, or tools.

1.1.4 Speed, temperature, and time are controlled with a push of a button and are displayed on the LED panel.

1.1.5 Over temperature and over speed controls ensure safety, long life, and sample integrity.

1.1.6 Advanced motor isolation system minimizes noise and vibration.

1.2 Design: Incu- Shaker Mini

1.3 Performance Specifications:

Specifications	E11920
Voltage	220-240V
Frequency	50-60 Hz
Temp. range	Ambient +5 to 70°C
Platform Dimensions	235x290mm
Speed	30 to 300 rpm
Overall dimensions	28x41x33cm



Shiv Ganga Industrial Estate,
Plot no.1 Lakeshwari, Bhagwanpur,
Roorkee – 247667 Uttarakhand, India

An ISO 9001 ,13485, 14001 & CE Certified

A Division of Abdos

Regd. Office:
13A, Government Place (East)
Kolkata – 700069, India

GST : 05AAGCA8605A1Z7

P : +91 (0) 1332 235277

E : labtech@abdosindia.com

W: www.abdoslifesciences.com

CIN : U85100WB2008PTC123903

Technical Data Sheet

Shaking Orbit	19mm
Capacity	20L
Weight	11.5kg

1.4 Material: Body – ABS

1.5 Color: Grey

1.6 E11920 – Incu-shaker mini with non-slip rubber mat - 9.5x11.5 inch

1.7 Spare Accessories for E11260, E11290, E11920

Specifications	
E11925	Universal magnetic platform for flask and tube racks 9.5x11.5 in
E11928	Universal Magnetic Clamp for 125ml , 250ml & 500ml flask
E11929	Magnetic Clamp for 1000ml
E11930	Magnetic Clamp for 500ml
E11931	Magnetic Clamp for 250ml
E11932	Magnetic Clamp for 125ml
E11933	Magnetic Clamp for 50ml
E11934	Magnetic Clamp for one Microplate
E11935	Magnetic Clamp for tube rack 30x15ml & 20x50ml tubes
E11936	Magnetic Clamp for tube rack 72x14mm culture tubes
E11937	Magnetic Clamp for tube rack 32x15ml tubes
E11938	Magnetic Clamp for tube rack 15x50ml tubes
E11939	Magnetic Clamp for tube rack 4x500ml or 600ml conical tubes
E11940	Magnetic Clamp for tube rack 12x1.5/2.0ml tubes horizontal
E11941	Magnetic Clamp for tube rack 4x15ml tubes horizontal
E11942	Magnetic Clamp for tube rack 2x50ml tubes horizontal
E11943	Magnetic Clamp for tilted holder for 15ml and 50ml tube rack
E11944	Universal Spring Platform



Shiv Ganga Industrial Estate,
Plot no.1 Lakeshwari, Bhagwanpur,
Roorkee – 247667 Uttarakhand, India

An ISO 9001 ,13485, 14001 & CE Certified

A Division of Abdos

Regd. Office:
13A, Government Place (East)
Kolkata – 700069, India

GST : 05AAGCA8605A1Z7

P : +91 (0) 1332 235277

E : labtech@abdosindia.com

W: www.abdoslifesciences.com

CIN : U85100WB2008PTC123903

Technical Data Sheet

E11948	Platform for 5x500ml
E11949	Platform for 8x250ml
E11950	Platform for 1x125ml
E11951	Platform for 4 microplates


Spare Accessories -

E11925 - Universal Magnetic Platform for Flasks and Tube racks (For E11920) - 9.5x11.5 inch

1.8 CE certified

1.9 Chemical resistance: see website (www.abdoslifesciences.com)

Warranty: 12 months from the date of Installation or 15 months from the date of Invoicing.

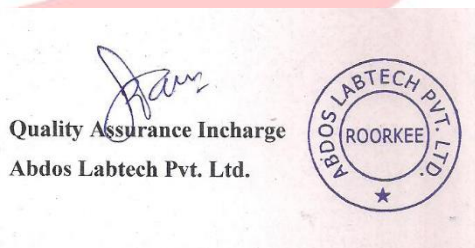

Shiv Ganga Industrial Estate,
Plot no.1 Lakeshwari, Bhagwanpur,
Roorkee – 247667 Uttrakhand, India

An ISO 9001 ,13485, 14001 & CE Certified

A Division of Abdos

Regd. Office:
13A, Government Place (East)
Kolkata – 700069, India

GST : 05AAGCA8605A1Z7

P : +91 (0) 1332 235277

E : labtech@abdosindia.com

W: www.abdoslifesciences.com

CIN : U85100WB2008PTC123903